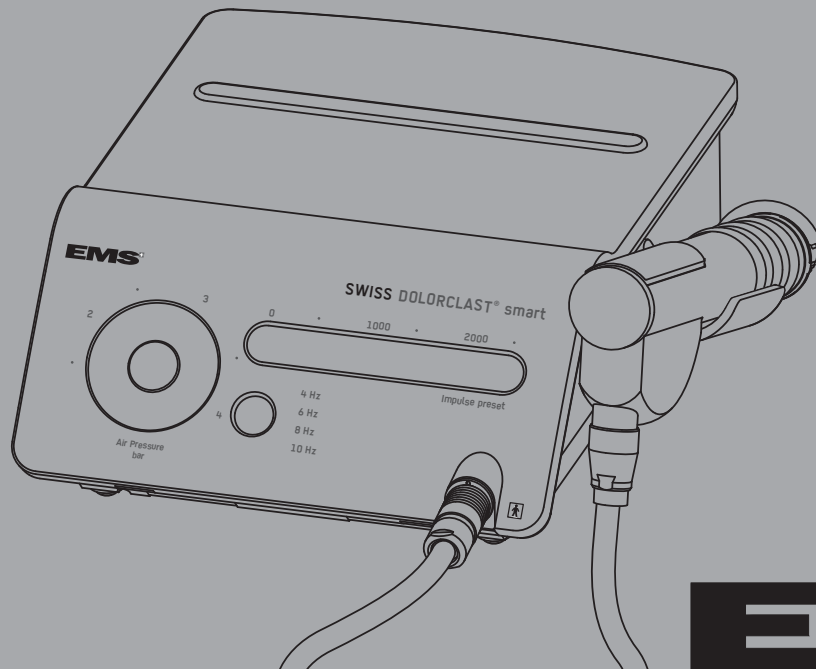
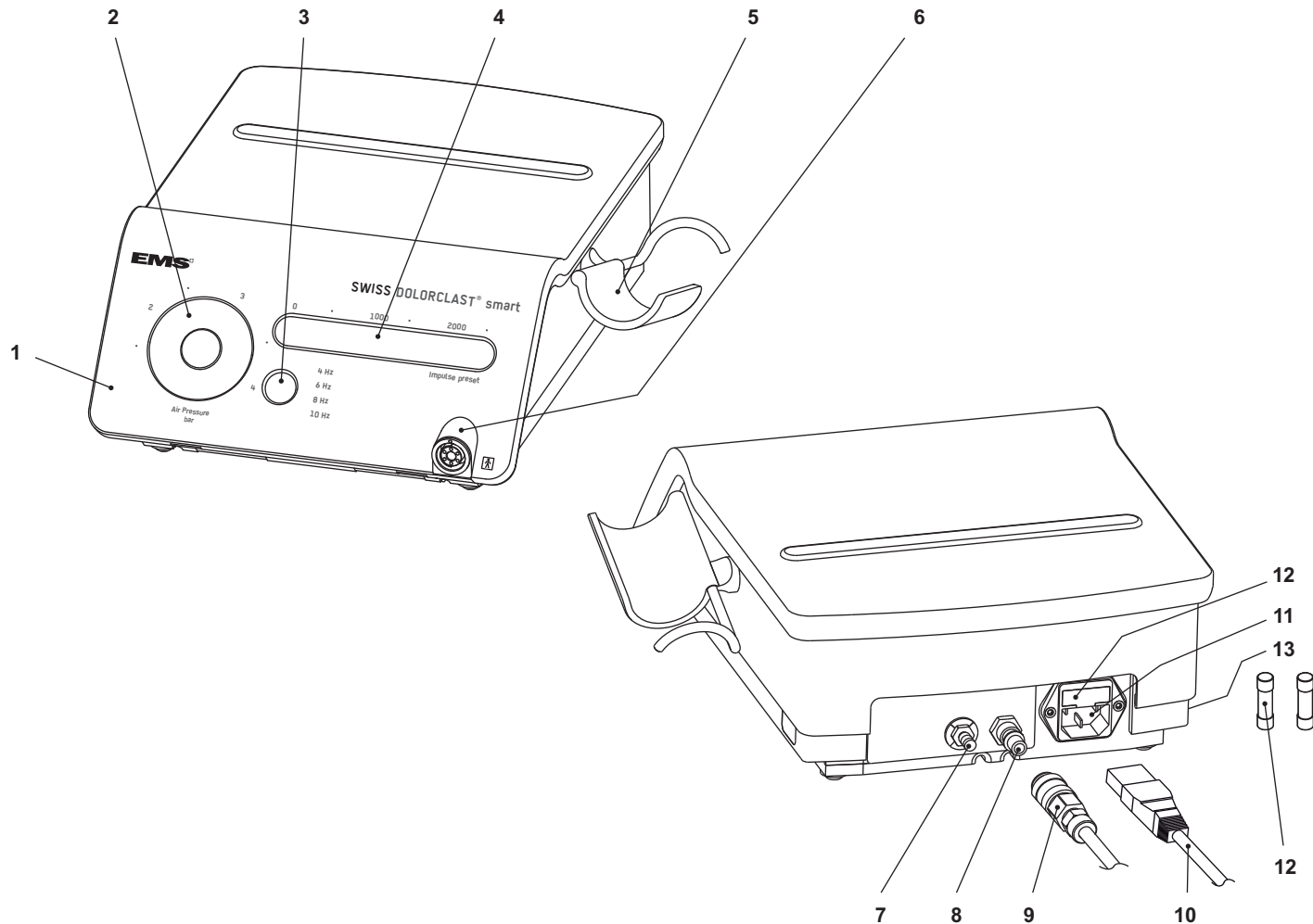
OPERATION INSTRUCTIONS
MODE D'EMPLOI
BEDIENUNGSANLEITUNG
INSTRUCCIONES DE EMPLEO
ISTRUZIONI PER L'USO
BRUKSANVISNING

SWISS DOLORCLAST[®] smart

Rx ONLY Caution! Federal (USA) law restricts this device to sale by or on the order of a physician

EMS 
ELECTRO MEDICAL SYSTEMS

COMPONENTS

- 1 Unit
- 2 Adjustment of working pressure
- 3 Frequency switch
- 4 Adjustment of number of pulses
- 5 Handpiece holder
- 6 Connector for the handpiece
- 7 Equipotential connector
- 8 Air hose connector
- 9 Compressed air tube
- 10 Mains power cord
- 11 Power socket
- 12 Fuses
- 13 Power switch

EMS supplies units with various accessories. The "Packing list" shows exactly what your unit includes.

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DEAR CUSTOMERS,

Thank you for buying a new EMS product. It meets the highest quality and safety standards.

The following indications have been successfully treated with the EMS Swiss DolorClast:

- Calcifying tendinitis of the shoulder
- Subacromial pain syndrome
- Tennis elbow (epicondylitis humeri radialis)
- Greater trochanteric pain syndrome
- Patella tip syndrome
- Rugby knee (Osgood-Schlatters disease)
- Medial tibial stress syndrome
- Insertional Achilles tendinopathy
- Midbody Achilles tendinopathy
- Plantar fasciitis (or fasciopathy)
- Myofascial pain syndrome (trigger point therapy), idiopathic low back pain and pseudoradicular syndromes
- Peyronie's disease (induratio penis plastica)
- Acute and chronic soft-tissue wounds (including diabetic foot ulcers)
- Primary and secondary lymphedema
- Cellulite

Several scientific publications and empirical reports have proven clinical efficacy and safety of the EMS Swiss DolorClast method for these indications.

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Contraindications:

- Treatment over air-filled tissue (lung, gut)
- Treatment of pre-ruptured tendons
- Treatment of pregnant women
- Treatment of patients under the age of 18 (except for the treatment of Osgood-Schlatters disease)
- Treatment of patients with blood-clotting disorders (including local thrombosis)
- Treatment of patients treated with oral anticoagulants
- Treatment of tissue with local tumors or local bacterial and/or viral infections
- Treatment of patients treated with cortisone

Risks:

- Pain and discomfort during and after treatment (anesthesia is not necessary)
- Reddening of the skin
- Petechia
- Swelling and numbness of the skin over the treatment area
- Potential damage of lung and gut tissue
- Potential rupture of pre-ruptured tendons
- Potential damage of immature epiphyses
- Possible dissemination of malignant tumors
- Potential damage of articular cartilage

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ABOUT THIS MANUAL

Note that the English version of this manual is the master from which translations derive. In case of any discrepancy, the binding version is the English text.

These operating instructions are to ensure the correct installation and use of this product.

Please read these operating instructions carefully as they explain all the most important details and procedures. Please pay special attention to the safety precautions.

Always keep these instructions close at hand.

To prevent injury to people and damage to property, please take note of the relevant warnings and remarks. They are marked as follows:

Danger
Risk of injury

Caution
Risk of property damage or environmental harm

Please note
Useful additional information and hints

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Förbjudet

Tillåtet

We would be pleased to answer your questions or listen to your suggestions. We do, of course, provide support in case of technical problems. Please call your approved EMS repair center or contact your dealer directly.

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EMS

CONTENTS

Installation	2
Connecting to the mains power	2
Connecting to the compressed air supply.....	4
Operating the unit	6
Connecting the handpiece	6
Switching On.....	8
Selecting the frequency.....	10
Setting the pressure	12
Presetting the number of pulses	14
Using the Swiss DolorClast smart	16
Cleaning and disinfecting.....	18
Cleaning the unit.....	18
Cleaning the handpiece.....	20
Maintenance.....	22
Fuses	22
Error mode.....	24
Display of software version	26
Safety precautions	28
Storing the product	30
Product disposal	32
Warranty.....	32
Accessories.....	32
EMS-Service.....	34
Symbols	36
Technical data	44
Electromagnetic compatibility.....	48
Troubleshooting	50

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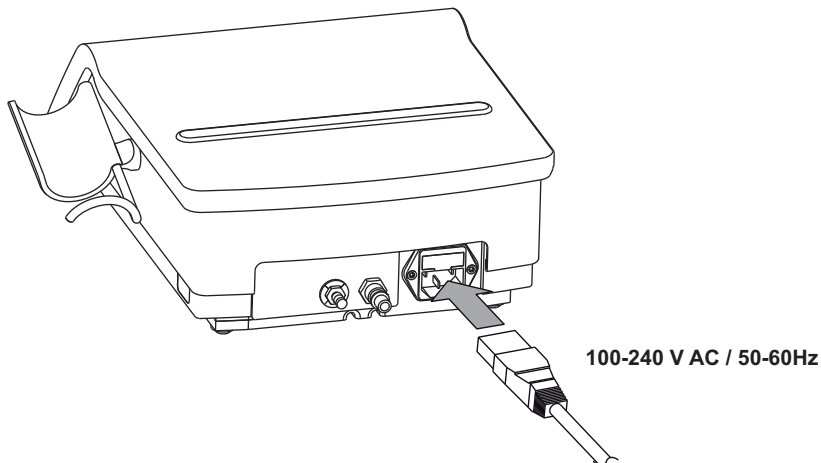
INSTALLATION

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INSTALLATION

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 To avoid risk of electric shock, this equipment must only be connected to a supply mains with protective earth.

 To prevent damage to the unit, check that its rated voltage corresponds to the local line voltage.

Use the unit only in a FI protected mains supply (FI = Residual current protection).

 Only for USA or Canada: Connect only to Hospital Grade receptacle.

 When applicable connect the equipotential plug to other medical devices used.

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SAFETY PRECAUTIONS

EMS and the distributor of this product accept no liability for direct or consequential injury or damage resulting from improper use, arising in particular through non-observance of the operating instructions, or improper preparation and maintenance.

 For further details on treatments, please refer to the operating instructions "Medical applications" for this product.

 Use for the intended purpose only. Before using the product, make sure that the operating instructions have been studied. This also applies to any equipment used with this product. Failure to strictly follow the operating instructions may result in serious injury to the patient or user or the product being damaged, possibly beyond repair.

 Extended use should be avoided in order to prevent the harmful effect of vibrations on the handpiece.

 Risk of explosion: do not use this product in the presence of inflammable anaesthetics or gases.

 This device must only be used by or on the order of a physician.

 In cases of tinnitus, it is recommended that the patient (treatment on the upper body) and the user (extended use) wear ear protection.

 Always examine the product for damage before commencing treatment. Damaged accessories or a damaged unit must not be used and must be replaced. Use original EMS spare parts and accessories only.

 Repairs should only be carried out by an approved EMS repair center.

 This product has been investigated with regard to safety from electrical shock and fire hazard. CSA (Canadian Standard Association) has not investigated other physiological effects.

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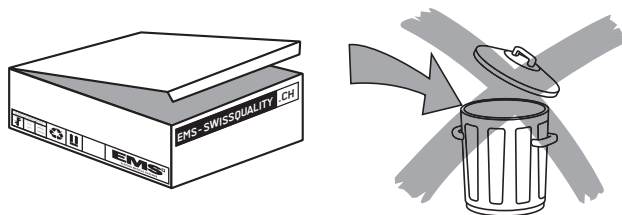
STORING THE PRODUCT

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 Keep the original packaging until the product is to be disposed of permanently. You can use it for shipping and storage.

Should you wish to store your product for an extended period:

- Proceed as described in the chapter "Cleaning and disinfecting".
- Place the product and all accessories in the original packaging.

Please refer to the storage and transport conditions in the "Technical data" section.

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PRODUCT DISPOSAL

The product must not be disposed of as domestic waste.

Should you wish to definitively dispose of the product, please comply with the regulations that apply in your country.

Waste Electrical and Electronic Equipment (WEEE) belonging to customers located in the European Union may be shipped to EMS for recycling in accordance to the WEEE regulations. The costs of recycling, exclusive of shipping fees, are covered by EMS.

WARRANTY

Damages due to non-adherence to the operating instructions or wear of parts are excluded from warranty.

The warranty of your product will be void if you try to open it.

ACCESSORIES

Accessories are available from EMS or an authorized dealer. Please contact your EMS customer service representative directly.

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EMS-SERVICE

Should your product need additional servicing or repairs, please send it to your dealer or directly to your approved EMS repair center.

In case of non-authorized repairs or damages due to non-adherence to the operating instructions, EMS accepts no liability whatsoever. This will also void the warranty.

No modification of this equipment is allowed.

Please clean, disinfect and sterilize the product and all accessories as described in the operating instructions before dispatching them.

It is best to ship your product in the original packaging. It protects your product against damage during shipment.

When sending the product directly to the approved EMS repair center, please include the name of your dealer. This simplifies processing.

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CSA marking with "C" identification for products in conformance with Canadian standards and "US" for products in conformance with US standards

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Disposal of Old Electrical & Electronic Equipment (Applicable in the European Union and other European countries with separate collection systems)

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TECHNICAL DATA

Description

Manufacturer	EMS SA, CH-1260 Nyon, Switzerland
Model	SWISS DOLORCLAST® smart
Classification EN 60601-1	Class I IP 20, unit Mode: Continuous operation
Applied part, Type BF	Swiss DolorClast handpiece
Classification 93 / 42 EEC	Class IIb
Supply voltage	100 V AC / 240 V AC
Mains frequency	50 - 60 Hz
Power consumption	40 VA
Fuses	Ø 5 x 20 mm, 250 V AC, time lag (Phase and Neutral) T 2.5 A L
Compressed air supply	5 - 6 bar (5000 - 6000 hPa) 20 ml/min at 4 bar
Operating mode	Multiple pulse mode repetition rate: 4 - 6 - 8 - 10 Hz
Operating pressure	1.5 - 4 bar
Weight	Approx. 2.9 kg
Dimensions in mm (HxWxD)	108 x 280 x 300
Operating conditions	+10°C to +40°C 30% to 75% relative humidity
Storage and transport conditions	-10°C to +40°C 10% to 95% relative humidity 500 hPa to 1060 hPa air pressure

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ELECTROMAGNETIC COMPATIBILITY

 This product needs special precautions regarding electromagnetic compatibility and needs to be installed and put into service according to the note "Electromagnetic compatibility".

 Certain types of mobile telecommunication equipment could potentially interfere with this product. The separation distances recommended in the note "Electromagnetic compatibility" must be taken into account.

 This product should not be used adjacent to or stacked with another unit. If adjacent or stacked use is necessary, this product should be observed to verify normal operation in the configuration in which it will be used.

 The use of accessories and cables other than those specified or sold by EMS as replacement parts, may result in increased emissions or decreased immunity of this product.

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TROUBLESHOOTING

The DolorClast smart provides an error code system in case of malfunction or incorrect handling:

1 LED blinking = Handpiece error; 2 LEDs blinking = Air pressure error; 3 LEDs blinking = Software error

Symptom	Error code	Solutions
The display does not light up when the unit is switched on	--	• Check the different electrical connections • Check the mains electricity supply • Check/replace the fuses and replace them with the type indicated on the back of the unit • Send the unit to your authorized EMS repair center
The display is blinking when the unit is switched on	--	• Send the unit to your authorized EMS repair center
Handpiece error	1 LED	• Check the different connections of the unit • Check the connection of the handpiece • Test the unit with another handpiece
Air pressure error	2 LEDs	• Check the different air connections • Check the air pressure supply is within the values indicated on the back of the unit • Check the compressor • Check the connection between the handpiece and unit • Test the unit with another handpiece • Send the unit to your authorized EMS repair center
Pressure not adjustable	--	• Switch OFF/ON the unit • Check the different air connections • Check the air pressure supply is within the values indicated on the back of the unit • Check the compressor • Send the unit to your authorized EMS repair center
Pressure reduction during operation	Repetitive "beep"	• Check the different air connections • Check the compressed air supply • Contact your authorized EMS repair center
Software error	3 LEDs	• Send the unit to your authorized EMS repair center

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EMS Corporation
EMS Electro Medical Systems GmbH